



**3-in-1:**

- Dissector
- Retractor
- Sizer



**The Triple Tool**

*Designed and Inspired by*  
**Dr. Luis Da Cruz**

For use in Gluteal Augmentation utilizing intra-muscular and subfascial techniques


Rotary dissection creates a safe, thick pocket in less than two minutes

Proprietary design facilitates control and precision

**Table for Implant Sizes**


Volume of Implant (cc)	Length of "B" (cm)*
280	16
320	17
430	18

\*Dimensions will change depending on the style and manufacturer.

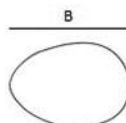

# MARINAMEDICAL

The New Wave in Gluteal Augmentation

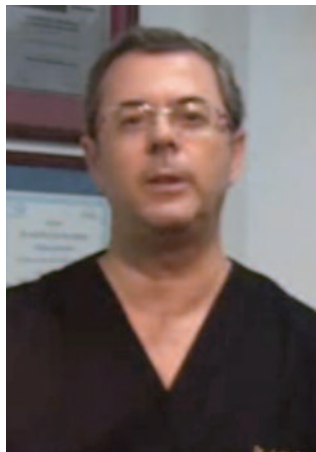

Dr. Luis Da Cruz

"Marina Medical is keen on innovation and is well-known as a pioneer in the plastic surgery field. I am **Dr. Luis Da Cruz**, and together we have developed a new device designed for gluteal augmentation with intramuscular and subfascial implants. Our goal is to standardize the steps in order to consistently achieve the appropriate depth for implant location and the correct thickness of gluteal muscle covering the implant.

The Triple Tool we developed is a unique 3-in1 instrument. It is a:  
**• Dissector • Retractor and • Sizer**  
**all in one instrument."**

**Marina Medical** is proud to offer superior instrumentation designed and inspired by Dr. Da Cruz.




*Inverting the position of the blades allows you to adapt to both gluteal sides for right or left dissection.*



*Together, the blades are connected by a stable handle and a pivot that creates a very strong union and helps make the displacement of tissue and opening of the full pocket for revision as a retractor.*



*The ergonomic design accommodates all anatomical implants.*

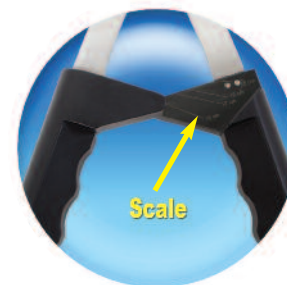

*A scale is located on the distal end of the handle.*



*Separately, each blade becomes a dissector and a retractor.*



*The thin blades with rounded edges allow gentle dissection of muscle fibers. The weight of the instrument also helps to provide better control.*



*The weight and every curve of the instrument plays a significant part in the stability and control provided.*



**Watch Dr. Da Cruz's video on the Marina Medical website for more information on the Triple Tool.**